



Coblation®: Better Patient Outcomes with Minimal Pain and Fast Recovery



New Coblator II


Introducing the Coblator® II Surgery System

- Variable Coagulation for *Improved Hemostasis*
- Lighter Patient Cable for reduced surgeon strain
- New advanced Foot Pedal for greater user efficiency
- Built-in Flow Control for set-up convenience

The ArthroCare® ENT Coblator II surgery system brings a new level of control and functionality to otolaryngology. Using patented Coblation technology, the ArthroCare ENT system removes tissue via molecular disintegration. For tissue coagulation, the new adjustable coagulation with its bipolar delivery of thermal energy provides localized therapy quickly and effectively.

Coblator II Surgery System


- Ⓐ *Controller*
- Ⓑ *Reusable Patient Cable*
- Ⓒ *Flow Control Valve*
- Ⓓ *Foot Control*
- *User's Manual*

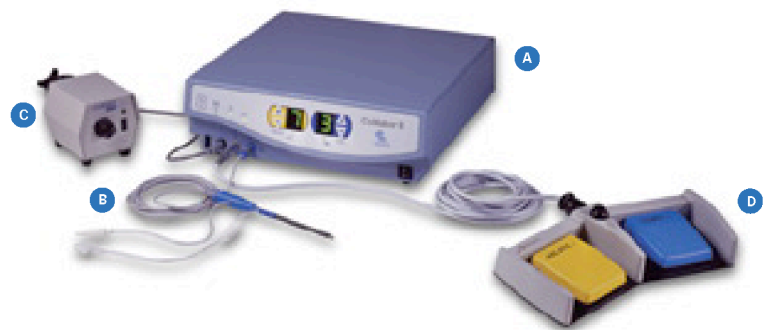

Coblator® II Family of Plasma Wands™

Available in various configurations, there is a Plasma Wand perfectly suited for several of the most commonly performed ENT surgical procedures. Wand options include built-in suction for soft tissue removal and enhanced visibility, tissue dissection and submucosal tissue channeling.

Suction Wands

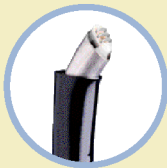

Suction Wands simultaneously dissect, ablate and remove tissue. The central suction channel removes debris and blood for a clear surgical view. The EVac™ 70 provides built-in suction and saline delivery.

Applications:

- Coblation® Assisted Tonsillectomy
- Adenoidectomy

Tissue Excision Wands



Tissue Excision Wands provide precise dissection using bipolar technology for minimal damage to surrounding tissue. The Plasma Hook includes built-in suction and saline delivery.

Applications:

- Uvulopalatoplasty (*RAUP*)

Channeling Wands



Channeling Wands combine controlled tissue removal and effective coagulative lesion formation. For submucosal tissue shrinkage in the nasal cavity, select the ReFlex Ultra™ 45. For submucosal tissue shrinkage in the oral cavity, or to treat snoring, choose the ReFlex Ultra™ 55.

Applications:

- Reduction of inferior turbinate (*ReFlex Ultra 45*)
- Soft Palate for the treatment of snoring (*ReFlex Ultra 55*)

System & Components	ITEM	CATALOG NO.	Coblator II Wands	ITEM	CATALOG NO.
	Coblator II	EC8000-01 <i>(Patient Cable not included)</i>		Suction EVac 70 Plasma Hook	E5870-01 E2850-01
Flow Control Valve	EA8002-00	Channelling ReFlex Ultra 45 ReFlex Ultra 55 Hummingbird®	E4845-01 E4855-01 E3845-01		
Patient Cable	EA8000-00				

IMPORTANT NOTE: Coblator II Wands can only be used in conjunction with the Coblator II controller and are not compatible with the original Coblator system.



ArthroCare ENT
680 Vaqueros Avenue
Sunnyvale, CA 94085-3523
toll-free **800-797-6520**
phone 408-736-0224
fax 408-736-0226
www.arthrocareENT.com

ArthroCare Europe AB
Baggensgatan 25
111 31 Stockholm
Sweden
phone +46 8 546 172 00
fax +46 8 546 172 39
info@arthrocare.se

CAUTION: Federal (USA) law restricts this device to sale by or on the order of a physician.

© 2004 ArthroCare Corporation. ArthroCare, Coblation, EVac, Coblator, ReFlex Ultra, and Hummingbird are trademarks and/or registered trademarks of ArthroCare Corporation.